


**INFORMATIVE NOTE**
**Name: REFLEXX 98**
**Date: 20/03/2024 Rev: 011**
**MANUFACTURER**

 REFLEXX SpA  
 Via Passeri, 2 46019 VIADANA (MN)  
 Tel.: + 39 0375 758891 Fax + 39 0375 464504  
[info@reflexx.com](mailto:info@reflexx.com)  
[www.reflexx.com](http://www.reflexx.com)

**Imballo - Pack**

50		10	
Taglia Size	Articolo Article	Taglia Size	Cartone Carton
R98 S	8032891631836	S	8032891636831
R98 M	8032891631843	M	8032891636848
R98 L	8032891631850	L	8032891636855
R98 XL	8032891631867	XL	8032891636862

**DESCRIPTION:** ■ disposable examination glove powder free, not sterile, high risk ■ **Material:** natural latex ■ **Shape:** ambidextrous; five fingers, beaded cuff, long cuff ■ **Colour:** blue ■ **Feature:** internally chlorinated; textured surface and textured finger ■ **Packing:** 50 pieces/box; 10 boxes/carton ■ **ISO 2859:** AQL 1.5 (Level G1) freedom from holes; AQL 2.5 (level G1) major defects; AQL 4 (level G1) minor defects

**CONFORMITY AND PERFORMANCES**
**Personal Protective Equipment cat III**

0465 ANCI\_CIMAC, Via Aguzzafame, 60/b 27029 Vigevano (PV)

Reg EU 2016/425 - EN ISO 374-1:2016 TYPE A - EN ISO 374-5:2016 - EN 374-2:2014 - EN 420:2003 + A1 2009

**Medical Device**

Reg EU 2017/745 - EN 455 1,2,3 &amp; 4

(CND T010201 NON SURGICAL LATEX GLOVES)



Reg CE 1935/2004 – Reg UE 10/2011 - DM 21/03/73. To such use carefully read the declaration of suitability for food contact (DICA).

DIMENSIONS		
Cod./Size	Length. min(mm)	Width. (±5 mm)
R98/ S (7)	280	85
R98/ M (8)	280	95
R98/ L (9)	280	105
R98/ XL (10)	280	115
R98/ XXL (11)	280	125

Thickness mm	Palm (+/- 0.06)	0.33
--------------	-----------------	------

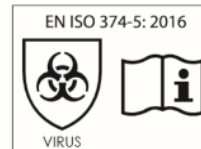
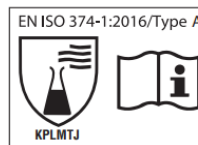
EN 455-2:2009/A1	Tensile strength	≥6 N
------------------	------------------	------

EN ISO 374-5:2016	Protection from microorganisms (bacteria and fungi)	Virus protection ISO 16604:2004
	pass	pass

**EN ISO 374-1:2016**

Chemical compound:	COD	EN 16523-1:2015 TYPE A		EN 374-4:2013
		Performance level	Breakthrough time	% degradation
Sodium hydroxide 40%	K	3	>60'	13,2%
Hydrogen peroxide 30%	P	3	>60'	29,7%
Sulfuric acid 96%	L	2	>30'	50,4%
Nitric acid 65%	M	2	>30'	36%
Formaldehyde 37%	T	3	>60'	17,9%
n-Heptane	J	2	>30'	64,7%

REG EU 2016/425


**EN 420:2003 + A1:2009**
**RISKS FROM WHICH THE DPI IS DESIGNED TO PROTECT:**

COD		Protection
K	Sodium hydroxide 40% CAS:1310-73-2	Protection from use of inorganic bases
P	Hydrogen peroxide 30% CAS: 7722-84-1	Protection from use of peroxides
L	Sulfuric acid 96% CAS: 7664-93-9	Protection from use of inorganic mineral acids, oxidants
M	Nitric acid 65% CAS: 7697-37-2	Protection from use of mineral acids
T	Formaldehyde 37% CAS: 50-00-0	Protection from the use of aldehydes
J	n-Heptane CAS: 142-82-5	Protection from the use of saturated hydrocarbons

**RECOMENDATIONS OF USE** Multi-purpose glove to protect the hand and/or the handled product, can also be used for general uses ■ The glove is suitable for contact with pure chemical products for a period of less than the permeation time indicated in the chemical resistance table ■ Results obtained under laboratory conditions may not reflect actual working time

**WARNINGS** Before using the product, try it in the real condition of use ■ The glove is a device subject to wear, so frequent replacement is recommended and in all cases where imperfections are observed ■ Wear gloves with dry, clean hands ■ The product contains natural latex which can cause allergic reactions of type I and IV hypersensitivity. In case of reaction, discontinue use and consult a doctor ■ The user must inform the manufacturer and the relative authority of any incidents occurring during the use of the protective equipment ■ Do not use gloves when there is a risk of getting caught in the moving parts of machinery

**STORAGE AND DISPOSAL** Store the gloves in their original packaging in a dry place. ■ Do not expose to direct sunlight and heat sources. ■ The expiry date is 3 years from the date of manufacture indicated on the sales package, and is valid for the product properly stored ■ Gloves must be disposed of according to local regulations ■ Declaration of Conformity on [www.reflexx.com](http://www.reflexx.com)